


## Enjoy the remarkable benefits of DermaNew

1/ Cleanse face remove any makeup, rinse your face and pat dry

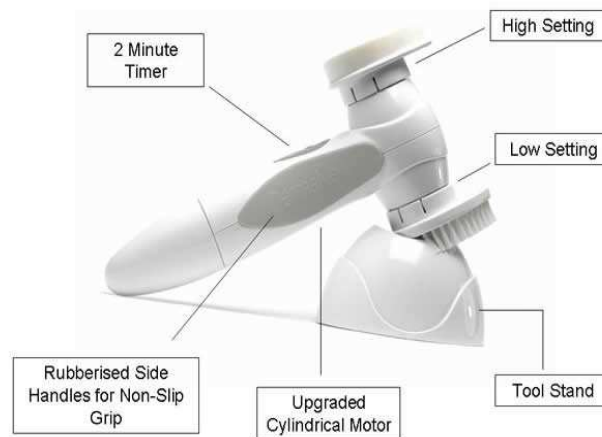

2/ Clip the Teardrop Foam Applicator onto the Dual Action Resurfacing Tool

One side is High setting for normal use and Low for use around delicate areas (eye area, neck, décolletage) or on sensitive skin.

It should be attached with the thinner teardrop end toward the top. This teardrop allows you to treat tighter areas of the face such as the bridge of the nose and around the eyes. (small foam applicator for the face & large applicator for the body)

3/ Apply a generous amount of the appropriate microdermabrasion cream according to you skin type or condition onto the foam head.

4/ Dot the Cream onto the areas to be treated, the forehead, cheeks, chin, nose and neck.



5/ Press the 'on' button once it will run for 2 minutes, by pressing it twice it will bypass the two minute auto shut of timer. To stop the device press it once more.

The dual action tool is waterproof, so can be used in the shower and also has a rubber grip handle, which makes it more secure to hold.

6/ Work the product into the skin for 2-5 minutes depending on your skin's needs and the area of the body you are treating. Using firm, upward circular movements, take care to avoid getting any cream in your eyes.

7/ Rinse the face well with warm water and pat dry.



8/ Apply either a moisturising cream or sunscreen onto your skin, depending on the time of day.

Remove the reusable Foam Applicator from the Firming Resurfacing Tool and clean it with hot water and soap. Allow sponge to dry before reusing. Replace batteries often for maximum results.

### **RECOMMENDED USE**

Use **DermaNew** 1-4 times a week as a regular addition to your skincare regime (morning or night)

DermaNew can be used on all areas of the face and body, including to treat callouses on feet, hands or elbows. (by using appropriate creams & attachments)

#### **Every day cleansing**

The DermaNew brush applicator can be used with the dual action tool and your desired facial cleansing product to give your skin a deep daily cleanse.

#### **CAUTION:**

Avoid eye area while cleansing or treating to avoid getting cleansers or lotions in your eyes. If eye contact occurs, flush gently and thoroughly with water. If irritation persists, consult physician. For external use only. If irritation occurs in other areas, discontinue use.