


Dermaroller®

Home-Care

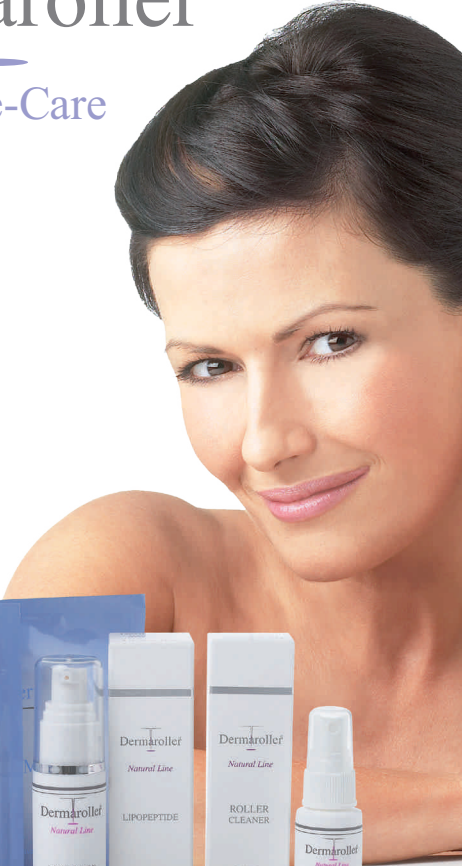

SMOOTHER, *BRIGHTER*,
YOUNGER
HEALTHIER LOOKING SKIN

Genuine Dermaroller™ Home-Care
Epidermal Skin Rejuvenation

GENUINE DERMAROLLER™ HOME-CARE EPIDERMAL SKIN REJUVENATION

The **Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation** systems can help improve the appearance of:

- Sun damaged & ageing skin
- Facial and décolleté fine lines and wrinkles

A specifically designed **Genuine Dermaroller™ HC9** device is used at home to stimulate the epidermis and to enhance the penetration of specific skincare products. The HC9 promotes smoother, brighter, younger, healthier looking skin.

**Only Genuine Dermaroller™ Devices
Incorporate Medizinisch-Nadel-Technik™
Technology and Genuine German
Medical Engineering**

Do not accept substitutes!

Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation

is only available from authorised clinics trained in

Genuine Dermaroller™ Therapy™



Quality, safety and efficacy assurance guaranteed.

WHAT IS GENUINE DERMAROLLER™ HOME-CARE EPIDERMAL SKIN REJUVENATION

Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation

uses a specifically designed home-care HC9 device that uses the same Medizinisch-Nadel-Technik™ medical-needle technology practitioners use with Genuine Dermaroller™ medical-care MC devices. The needles however are much shorter so they stay within the epidermal surface layer of the skin and the device is designed for usage, cleaning, storage and re-use safely and comfortably within your own home. Used on its own or as part of a system with Dermaroller™ dermo-cosmetic products, as described later in this leaflet, it will encourage smoother, brighter, younger, healthier looking skin.

Genuine Dermaroller™ Home-Care Epidermal Skin


is less invasive than and not a replacement for the **Genuine Dermaroller™ Therapy Skin Rejuvenation and Scar Repair** medical-care procedures carried out in the clinic*. It is however a perfect complement to help enhance and prolong results achieved with the in-clinic procedures or you may agree with your practitioner that **Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation** without the need for medical-care in-clinic procedures is currently the best approach, the latter is most-likely to be the case if you have younger or less damaged skin.

How it Works

The **Genuine Dermaroller™ HC9** when used as described in the 'How to Use' section creates thousands of tiny "medizinisch-nadel" needle columns in the epidermis. This physical stimulus triggers the epidermal skin cells to release growth factors which increase epidermal cell proliferation and dermal renewal. It has also been shown to stimulate collagen production and proteins that help improve the structure of the skin and as a result help it look smoother, brighter and healthier. The thousands of tiny "medizinisch-nadel" needle columns remain open for several minutes after a procedure, these can help aid the penetration and action of anti-ageing and skin enhancing dermo-cosmetic products. These columns will then close over enabling your stimulated epidermis to provide a healthy barrier.

*To find out more about **Genuine Dermaroller™ Therapy™**
visit: www.genuinedermaroller.co.uk

How to Use

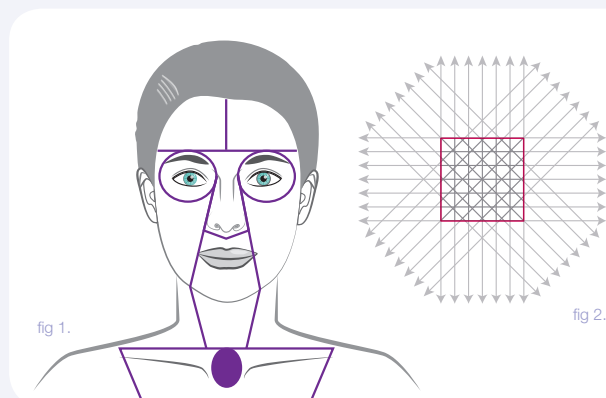

The **Genuine Dermaroller™ HC9** should be used initially twice a week on an evening. Firstly consider the face as sub-areas as demonstrated in figure 1 and make sure your skin is clean.

Working on one sub-area at a time use a firm grip so as to control the HC9 and make a pass forwards and backwards. Complete 3 passes in the same location and then repeat this process moving across the sub-area vertically, horizontally and diagonally. Avoid sideways movement and downward pressure.

If after 2 to 3 weeks use your skin is tolerating the use of the **Genuine Dermaroller™ HC9** the frequency of use and number of passes described above can be increased.

After each use place the HC9 device under a running hot tap and then carefully shake dry taking care not to knock the needles. Spray the Dermaroller™ Roller Cleaner on the HC9 roller head and needles and place and store in its container. The HC9 device has a use life of approximately 9 - 12 months with proper care and with an average use of 2 to 3 times per week. The **Genuine Dermaroller™ HC9** device should be used under the guidance and based on the advice of your aesthetic skin professional. It should not be used on active skin conditions such as acne, eczema or psoriasis. If in doubt consult a professional. If the use of the **Genuine Dermaroller™ HC9** causes irritation or excessive redness stop its use immediately.

**NEVER SHARE A
GENUINE DERMAROLLER™ HC9
WITH ANY OTHER USER!**



Further information available on www.genuinedermaroller.co.uk

GENUINE DERMAROLLER™ SKIN REGENERATING HC EPIDERMAL SYSTEM

The **Genuine Dermaroller™ Lipopeptide (Vitamin A) Serum**, used in combination with the **Genuine Dermaroller™ HC9** helps to regenerate ageing and sun damaged skin.

The **Dermaroller™ Lipopeptide (Vitamin A) Serum**; now combines high concentrations of the powerful skin regenerating Vitamin A, an advanced lipopeptide to boost skin enhancement, with the protecting Vitamin E anti-oxidant and enhanced plant-derived skin moisturising agents to further protect and hydrate the skin.



How to Use

Apply a thin layer and massage gently into the skin two nights a week until the skin is acclimatised after which it can be used more regularly. You may observe some redness and feel some “skin tingle” when first using the product. This is normal as the skin becomes acclimatised to the Vitamin A.

Use of a high protection UVA & UVB sunscreen during the day is essential during the period that the **Dermaroller™ Lipopeptide (Vitamin A) Serum** is used. AesthetiCare® recommends **Heliocare® Advanced SPF50 Gels**. Failure to do so could result in the skin developing areas of hyper-pigmentation.

Initially use the **Dermaroller™ Lipopeptide (Vitamin A) Serum** and **Genuine Dermaroller™ HC9** device on different nights. Once you do not get skin redness or tingle after the **Dermaroller™ Lipopeptide (Vitamin A) Serum** use you can try applying after using your HC9 device. First roll a test area where you have been using the product and apply a thin layer of the **Dermaroller™ Lipopeptide (Vitamin A) Serum** and wait for 24 hours. If it is well tolerated you can use on larger areas where the lipopeptide has been previously used and well tolerated. If excessive redness or skin tingle occurs revert to using on different nights. Use of the **Dermaroller™ Lipopeptide (Vitamin A) Serum** and **Genuine Dermaroller™ HC9** device as described above will help provide smoother, brighter, younger, healthier looking skin.

GENUINE DERMAROLLER™ SKIN HYDRATING HC EPIDERMAL SYSTEM

The **Genuine Dermaroller™ Mask** used in combination with the **Genuine Dermaroller™ HC9** device hydrates and enhances skin that is prone to dryness or in need of a hydrating boost.

The **Dermaroller™ Mask** is designed to intensively hydrate and replenish skin. It is a soft silky mask imbued with a novel serum containing the intensive skin hydrator Hyaluronic acid, the essential Ceramide skin lipid, the protecting anti-oxidant Vitamin E and the skin regenerating peptide Matrixyl.



How to Use

Following the use of the **Genuine Dermaroller™ HC9** as described earlier apply the **Dermaroller™ Mask** as instructed below:

- With clean hands, open the sachet, take out the mask, unfold it and remove the blue protective layer
- Cover the face with the silky side down. Remove the white protection layer and pat the silk evenly on the skin
- Leave it on the skin for a minimum of 5 minutes; maximum time is when the silk starts to dry
- Remove the mask and gently massage the rest of the serum into the skin

The use of the **Dermaroller™ Mask** and **Genuine Dermaroller™ HC9** device as described above will leave the skin feeling and looking hydrated, refreshed and healthier. This can be repeated 2 to 3 times a week or as advised by your aesthetic skin professional.

If the use of the **Dermaroller™ Mask** leads to irritation or redness stop its use immediately and wash the affected area.

IT IS IMPORTANT TO FOLLOW THESE AND YOUR AESTHETIC
PROFESSIONAL'S GUIDELINES FOR USE.



Dermaroller®

Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation devices and products are only available from trained medical practitioners and authorised clinics. Ask your aesthetic professional about **Genuine Dermaroller™ Home-Care Epidermal Skin Rejuvenation** or visit www.genuinedermaroller.co.uk for further information.



FOLLOW US ON TWITTER

@dermarolleruk



FIND US ON FACEBOOK

www.facebook.com/genuinedermaroller

CLINIC DETAILS



AESTHETICARE®

The Professional Choice

Advanced evidence based dermatological
medical-devices, pharmaceuticals & medigrade skincare

freephone **0800 0195 322**
www.aestheticare.co.uk

Unit 124, Thorp Arch Estate, Wetherby,
West Yorkshire LS23 7BJ

A division of Ferndale Pharmaceuticals Ltd.