

PURE RAYZ

FULL-FACE PHILOSOPHY BY BABY QUASAR



INSTRUCTION MANUAL



PLEASE READ ALL INSTRUCTIONS BEFORE USE



TABLE OF CONTENTS

WHO WE ARE	1
HOW IT WORKS	1
LIGHTWAVE TECHNOLOGY	2
SETUP	3
LED LIGHT PATTERN	3
TEST FOR LIGHT SENSITIVE SKIN	3
TREATMENT REGIMEN	3
TREATMENT STEPS	4
TIPS FOR USE	4
SERUMS, CREAMS, LOTIONS	5
AFTER EACH TREATMENT	5
CARE AND MAINTENANCE	5
PROBLEMS	6
RISKS	6
WARNINGS	6
PRECAUTIONS	7
WARRANTY INFORMATION	7-8
PURE RAYZ TOOL DESCRIPTION	8
CUSTOMER SERVICE	9

WHO WE ARE / HOW IT WORKS

WHO WE ARE

*qua-sar (kway-sar): any class of celestial objects that resemble stars but whose large redshift and apparent brightness imply extreme energy output.

Quasar Bio-Tech was founded in 2005 on the belief that offering premium products with outstanding customer care was the only way to do business. Today, our continued success is due in large part to you—our clients. We would like to thank you for your continued support of our endeavors.

Located in Sarasota, FL since 2007, we pride ourselves in being a small business with big dreams to continue offering skin care products that aid our clients in feeling and looking their best. We are extremely proud to manufacture our products here in the United States.

Quasar Bio-Tech looks forward to your continued support, suggestions, and success stories.

HOW IT WORKS

FULL FACE PHILOSOPHY, SKIN THAT GLOWS FROM THE INSIDE OUT IN ONLY MINUTES A DAY.


Your body begins producing less collagen and elastin as we grow older. This leads to wrinkles and fine lines. Topical products can only reach so far, but the Pure Rayz light penetrates through the layers of skin to stimulate collagen and elastin development thereby diminishing fine lines and wrinkles.

Pure Rayz is a skin care tool that produces a gentle and natural light that works on all skin types and tones. Unique among skin care treatments, it produces long-lasting results and requires only occasional treatments to maintain benefits.

LIGHT WAVELENGTH TECHNOLOGY

LIGHT TECHNOLOGY

Pure Rayz combines four distinct wavelengths of light specifically designed to penetrate deeper than lotions and creams alone.

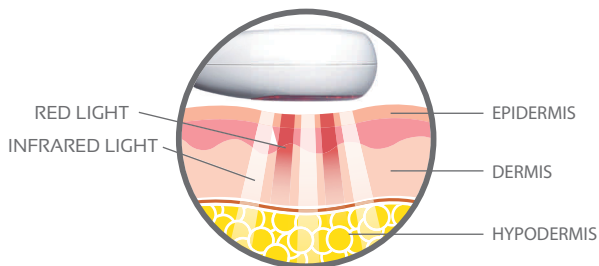

SETUP

1. Insert the power plug into the bottom of the Pure Rayz tool.
2. Plug the power supply into the wall.
3. Press the On/Off switch to make sure tool is functional.

PURE RAYZ LED LIGHT PATTERN

Note: The 10 white LEDs (5 on the right and 5 on the left) in the illustration are infrared LEDs and will appear to be off even though they are working. You cannot see infrared light with the naked eye.



TEST FOR LIGHT SENSITIVE SKIN

1. Place the Pure Rayz tool directly on your arm. Turn on the tool and move it slowly in the area for the full three-minute cycle.
2. If you have pain or redness, you may have light sensitive skin.

TREATMENT REGIMEN

Treatment is three (3) minutes in each area you are treating five (5) days per week for the first eight (8) weeks. The Pure Rayz tool is cleared for use up to five days per week for the first eight weeks. There is a risk, however, of drying and redness of the skin. If this occurs, please reduce treatments by one day per week until dryness subsides.

After the initial treatment period of eight weeks, you should treat once per week or as needed to ensure that collagen and elastin production remains active and that your skin retains its healthy, youthful glow.

TREATMENT STEPS / TIPS FOR USE

TREATMENT STEPS


1. Start by cleansing your face.
2. If you choose to use a serum, apply at this time.
3. Place the Pure Rayz tool gently on the skin in the first area you are treating.

Note: *Do not place the tool directly over the eyelids and keep your eyes closed when treating around the eyes.*

3. Move the tool in a slow circular or sweeping motion. Treat within the area for three (3) minutes until the tool emits a soft tone and the tool shuts off.
4. Turn the tool back on and repeat until all areas have been treated.

TIPS FOR USE

- When treating around the eyes, use a "C" motion around the left eye and a REVERSE "C" motion around the right eye.
- Place the tool on your eyebrow and move the tool in a slow "C" motion or reverse "C" motion around the eye (eyebrow, corner of eye and directly under the eye). This will treat the forehead, "11" lines, crow's feet, under the eye and the cheek.
- When treating the bottom of the face, use a left to right "sweeping" motion across the bottom of your face and directly over your lips.

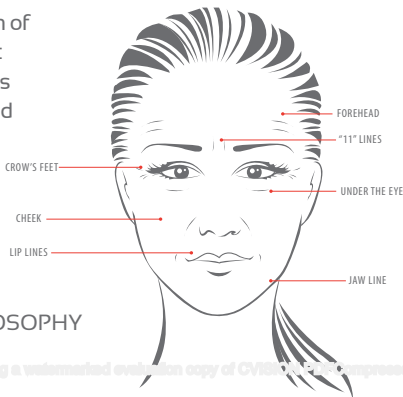

FULL FACE PHILOSOPHY

SERUMS, CREAMS AND LOTIONS

You can use the Quasar Anti-Aging serum or any serum before doing treatments with the Quasar tools. Serums are light and will emulsify quickly into the skin. We suggest you do a product light-sensitivity test on your arm before applying product to your face and doing treatments with the Quasar tools.

Refrain from using creams or lotions until AFTER you complete your treatment. Creams and lotions can cause a barrier to the skin preventing light penetration.

INJECTIONS/FILLERS/MICRODERMABRASION/CHEMICAL PEELS

After receiving injections/ fillers, wait 2 weeks before treating the area of the injection/filler. This will give the injection/filler time to settle.

You may use your Pure Rayz tool immediately after a microdermabrasion treatment or chemical peel.

NOTE: Make sure all products from the chemical peel are completely washed off the face before using the light.

AFTER EACH TREATMENT

We recommend using a moisturizer after each treatment to ensure your skin is sufficiently hydrated. You may apply make-up or other skin care products as you like. Quasar Tools do not cause sun sensitivity.

CARE AND MAINTENANCE

Use a soft, damp cloth or alcohol wipe to clean the tool's surfaces. Store the tool in a cool, dry place. Do not use any disinfectants (except an alcohol wipe) on your tool.

PROBLEMS / RISKS / WARNINGS

PROBLEMS?

***DO NOT TRY TO FIX THE PURE RAYZ TOOL YOURSELF.
CONTACT US IMMEDIATELY IF ANY PROBLEMS ARISE.**

RISKS

- There is a risk of light sensitivity.
- During treatments, there is a risk of drying and redness of the skin.
- Do not place the tool directly over the eyelids and keep your eyes closed when treating around the eyes.

Your safety is our primary concern. Read instruction manual thoroughly before using Pure Rayz tool.

WARNINGS

- DO NOT USE if Pure Rayz or Pure Rayz power supply has fallen in water or is damaged in any way to prevent risk of electric shock.
- DO NOT bend or break the wire from the power supply to prevent risk of electric shock.
- DO NOT USE any other power supply with your Pure Rayz, unless it came directly from Quasar Bio-Tech as a replacement cord, as it could cause a fire.
- DO NOT USE your Pure Rayz if you suspect it to be defective in any way and contact us immediately.
- DO NOT take apart your Pure Rayz for any reason to prevent risk of fire or damage to you or your tool.
- NOTE: Warranty is VOID if you, for any reason, disassemble your Pure Rayz.

PRECAUTIONS / WARRANTY INFORMATION

PRECAUTIONS

- ALWAYS unplug your Pure Rayz tool from the wall after use.
- DO NOT look directly into the light as it may cause pain or headache. If you know you have light-sensitive eyes, we recommend you wear protective eyewear when using the Pure Rayz tool.
- DO NOT USE Pure Rayz if your skin is sensitive to light, whether for genetic reasons or resulting from taking certain medications. If you are unsure of your skin's light-sensitivity, refer back to the TEST FOR LIGHT-SENSITIVE SKIN section.
- During treatment, do not use the Pure Rayz tool directly over the eyes or eyelids.
- DO NOT USE if you are pregnant or suspect you may be pregnant as no studies exist to confirm or deny possible damage. Always consult your physician for any health-related questions involving this tool and any specific conditions you may have.
- DO call a physician right away if you think you are experiencing a negative reaction of any kind after using the Pure Rayz tool, such as redness of the skin lasting longer than two hours after treatment.

WARRANTY INFORMATION

It is important to register your device with us to ensure the FIVE-YEAR warranty offered. It only takes a minute to register your tool at www.babyquasar.com or just complete the enclosed warranty card and drop in the mail. Registration must be completed and returned within 20 days of purchase. Your Pure Rayz tool is guaranteed free from defects in materials and workmanship under normal use

TOOL DESCRIPTION

for a period of FIVE (5) years after its original purchase date. This warranty is not transferable if you sell your Pure Rayz tool to another party. If any problems surface during this five-year period, contact our friendly, helpful customer service representatives immediately to allow us the opportunity to rectify any concerns.

North America, call toll-free (800) 944-1523 for repair or replacement.

Outside of North America, call (941) 306-5812.

Contact information also available at www.babyquasar.com.

PURE RAYZ TOOL DESCRIPTION

SPECIFICATIONS

- Size – 3.6 inches X 2.6 inches X 1.3 inches
- Weight – 2 ounces
- 3 minute treatment timer
- 100-240 vac, 50-60 Hz universal power supply



Quasar Bio-Tech, Inc. reserves the right to replace a defective product with the most comparable product currently available. Your exclusive remedy for any defective product covered by this five-year warranty is limited to the repair or replacement of the defective product. Quasar Bio-Tech's entire liability for any defective product shall not exceed the purchase price for the defective product. Manufactured proudly in the United States by: Quasar Bio-Tech, Inc. Sarasota, FL 34243 (941) 306-5812.

CUSTOMER SERVICE

CUSTOMER SERVICE

Our friendly and helpful customer service representatives are eager to help you in any way possible. Please feel free to call us from 9 - 5 p.m., Monday through Friday, Eastern Time, USA. You may also email us questions.

Please contact us:

Quasar Bio-Tech, Inc.
1431 Tallevast Rd. Sarasota, FL 34243
(800) 944-1523 (toll-free)
(941) 306-5812 (voice)
(941) 306-5816 (fax)
info@quasarbiotech.com
www.babyquasar.com

PURE RAYZ

FULL-FACE PHILOSOPHY BY BABY QUASAR

All Rights Reserved
Printed 05 September 2013

LAB-058-A